



SureSigns VM6

The versatility you need in a portable bedside monitor

PHILIPS

SureSigns VM6.

Your days can get pretty complicated. They are already full, and likely to become more complex as time goes on. Trust Philips to give you what you need to help care for your patients. The SureSigns VM6, for example, provides you with numerous essential measurements in one easy-to-use, compact monitor.

Every Philips SureSigns monitor blends vital measurement and monitoring capabilities with convenience and reliability. Patient monitoring—and quality patient care—have never been more clear.



Easy to learn and easy to love

The VM6, like all SureSigns monitors, is easy to set up, easy to learn, easy to use, easy to support, and, above all, easy to love. Measurements include user-selectable 3- or 5-lead ECG, pulse oximetry, non-invasive blood pressure, respiration, basic arrhythmia, temperature, and one invasive blood pressure (optional). Intuitive navigation puts your most commonly used features at your fingertips. The crisp color display is bright and easy to read, even in low-light situations. SureSigns VM6 is lightweight, yet durable, with a comfortable handle for easy mobility, and a battery life of 4 hours.

Connect to the SureSigns VSV: a simple way to be everywhere at once

Having a clear picture of each patient's status, at a single glance, can help relieve stress on care teams. The Philips SureSigns Vital Signs Viewer is an affordable, dependable, and simple way to add the benefits of central viewing and alarming without the complexity or expense of a traditional PC-based central station. Connect up to 12 SureSigns VM monitors at a time. It is one more way Philips is making it easier to free yourself to do exactly what you need to do to care for your patients.



All the support you need

The SureSigns VM6 is backed by high-quality technology and support. Software upgrades, via the USB port, are remarkably quick and easy. The sleek modular design simplifies set up and repair. Our warranty adds to your peace of mind. And all Philips monitors come with 24-hour telephone support from our global response centers*.

The Philips VSV allows you to check the status of up to 12 patients at a glance.

* Available only in U.S.

Beyond the basics.

Integrated bedrail hanger for maximum flexibility.

Large 8.4" (20.3 cm) color display with up to four waveforms and numerics is bright and easy to read.

Sturdy handle and rugged housing for easy portability.

Lightweight at only 3.2 kg (7.1 lbs).

Lithium ion battery is standard.

Fixed keys for quick access to main functions.

Non-invasive blood pressure (NBP), including manual, interval, and STAT modes.

Integrated strip recorder for convenient documentation (optional).

Up to 96 hours of measurement trends.

Bar code scanner for convenient entry of patient ID (optional).

USB port for easy software upgrades.

User-selectable 3- and 5-lead ECG with pacemaker detection, including a "leads off" indicator.

Philips SpO₂ technology compatible with a wide variety of adult, pediatric, and neonatal sensors from Philips, Nellcor,** and Masimo.**

** Adapter cable required.

| Product specifications | |
|------------------------|---|
| Width, height, depth | W: 26 cm (10.2"), H: 22 cm (8.6"), D: 14.5 cm (5.7") |
| Weight | 3.2 kg (7.1 lbs), excluding optional recorder |
| Screen | 8.4" (20.3 cm) SVGA TFT-AM LCD |
| NBP | Oscillometric using stepwise deflation pressure <ul style="list-style-type: none"> Adult/pediatric measurement range <ul style="list-style-type: none"> Systolic: 30 – 255 mmHg Diastolic: 15 – 220 mmHg Mean: 20 – 235 mmHg Neonatal measurement range <ul style="list-style-type: none"> Systolic: 30 – 135 mmHg Diastolic: 15 – 110 mmHg Mean: 20 – 125 mmHg Blood pressure accuracy: ± 5 mmHg Pulse rate accuracy: ± 2 bpm or $\pm 2\%$ (whichever is greater) |
| IBP (optional) | <ul style="list-style-type: none"> Range: -40 – 360 mmHg Accuracy: ± 4 bpm or $\pm 4\%$ (whichever is greater) |
| SpO ₂ | <ul style="list-style-type: none"> Meets EN ISO 9919 standard Measurement range: 0 – 100% Accuracy depends on sensor |
| ECG | <ul style="list-style-type: none"> 3- or 5-lead; heart rate range: 15 – 300 bpm Respiration rate range: 4 – 150 rpm |
| Temperature | <ul style="list-style-type: none"> Mode: continuous Method: esophageal, rectal, or skin Range: 25 – 45°C (77.0 – 113.0°F) Accuracy: ± 0.1°C |
| Battery | <ul style="list-style-type: none"> Lithium ion, 11.1V Operating time: 4 hours with continuous monitoring and NBP every 15 minutes |
| Data output | HL7 format, via Ethernet port |
| Patients | Adult, pediatric, neonatal |

| Supplies included* | To re-order |
|--|---|
| Adult antimicrobial NBP cuff | M4555A |
| NBP interconnect tubing | M1599B |
| Adult finger SpO ₂ sensor | M1196A |
| 3-lead ECG lead set with grabber cable | 989803143171 (IEC) or 989803143181 (AAMI) |
| Recorder paper (5 rolls)** | 989803136891 |

| Also available | To order |
|------------------------|--------------|
| Roll stand with basket | 989803144001 |
| Wall mount | 989803144011 |
| Bar code scanner kit | 989803147821 |

A simple way to be everywhere at once


Connect with the SureSigns Vital Signs Viewer to show status from up to 12 patient monitors at a time.

Getting started

Each monitor comes with set-up guide, instructions for use, quick reference cards, user training DVD, and service manual on CD.



Supplies for the SureSigns VM6 include a wide selection of disposable and reusable SpO₂ sensors, all designed for accuracy, performance, and durability.

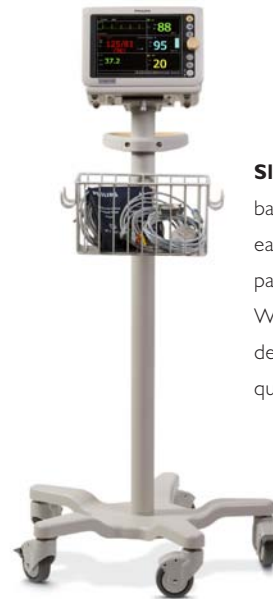

Sleek roll stand with integrated basket and locking casters makes it easier to take the SureSigns VM6 patient monitor where it's needed. Wall mount also available. Each is designed for easy mounting and quick release.

* Additional supplies available. Please contact your Philips representative for more information.

** Depending on model ordered.

| | |
|--------|------------------------|
| 863064 | SureSigns VM6 |
| 863065 | SureSigns VM6 with IBP |
| R01 | Recorder (optional) |
| 863067 | SureSigns VSV |

For more information, visit www.philips.com/SureSigns

buy online⁺
<http://shop.medical.philips.com>

+ Available only in U.S.



© 2008 Koninklijke Philips Electronics N.V.
 All rights are reserved. Philips Medical Systems Nederland B.V. reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

Nellcor is a registered trademark of Nellcor Puritan Bennett, Inc. Masimo is a registered trademark of Masimo Corporation.

Philips Healthcare is part of Royal Philips Electronics

www.philips.com/healthcare
 healthcare@philips.com
 fax: +31 40 27 64 887

Printed in The Netherlands
 4522 962 31551/862 * APR 2008

Philips Healthcare
 Global Information Center
 P.O. Box 1286
 5602 BG Eindhoven
 The Netherlands